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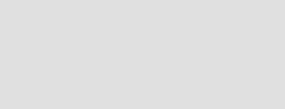
Instructions for Use

Indications

The XIVE® custom abutment ZrO₂ by Compartis® is a custom-made device fabricated **individually for each patient.**

The **XIVE® custom abutment ZrO₂ by Compartis®** has been approved for the fabrication of single-piece abutments in the anterior region serving to support highly aesthetic maxillary or mandibular single crowns at sites 13 to 23 and 43 to 43 (including cases with very low mucosal thickness) in connection with XIVE® CXz system implants. The **XIVE® custom abutment ZrO₂ by Compartis®** can be used in orthoaxial situations as well as in cases of custom implant angles up to 15°.

Due to the necessary level of precision of the implant interface geometry, we offer this indication only through our Compartis® integrated systems network production centres.

Important note: Implant-supported restorations are characterized by the absence of a tactile response on the part of the prepared tooth similar to Sharpey fibres as the masticatory forces are transmitted to the restoration. To prevent any damaged by peak loads, the abutments and frameworks as well as the ceramic veneer must be produced with meticulous care (uniform thickness of all ceramic layers and strict adhesion to the most recent firing recommendations for Cercon® ceram kiss and Cercon® ceram love).

Contraindications

- Use of the **XIVE® custom abutment ZrO₂ by Compartis®** where the length of the implant used is smaller than the length of the clinical crown to be replaced is contraindicated.
- Use of the **XIVE® custom abutment ZrO₂ by Compartis®** in patients exhibiting bruxism or function-resistant eror function is contraindicated.

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Warnings

Read these Instructions for Use completely before using the **XIVE® custom abutment ZrO₂ by Compartis®**. The **XIVE® custom abutment ZrO₂ by Compartis®** must be used only for the approved indications according to the general rules for dental/surgical procedures, taking all applicable occupational health and safety regulations into due account. The **XIVE® custom abutment ZrO₂ by Compartis®** has been approved for medical/dental use in conjunction with the XIVE® implant system only. The abutment must be used only in perfect condition. If you are uncertain about the indication or the correct mode of application, do not use the product until all open points have been resolved. We warrant the quality of our products within the framework of our Terms and Conditions for Sales and Deliveries. The following descriptions are not sufficient to ensure a professionally correct application of the abutment by treatment providers in experienced

in their use. The user must be thoroughly familiar with dental surgery and prosthodontics, including diagnosis and preoperative planning and/or laboratory procedures. Because the use of the product is beyond our control, we will not be responsible for any damage caused in conjunction with this. The responsibility rests entirely with the treatment provider.

- The following complications may be associated with the use of the **XIVE® custom abutment ZrO₂ by Compartis®**:
 - Aspiration or swallowing of components used in the patient’s mouth
 - Fracture of the screw in the abutment due to excessive torque
 - Rotation of the abutment due to insufficient torque
 - Fracture of the abutment due to overload or insufficient bulk following individual reduction
 - Damage to the abutment following reduction (even when irrigated)
 - Soft-tissue irritation or inflammation caused by cement residue

Precautions

Make sure to take the following precautions before and during each treatment session:

- Before each procedure, make sure that all components, instruments and auxiliary materials required are present in the necessary quantities, complete and functional.
- For your own safety, wear suitable protective clothing at all times.
- Position the patient so that the risk of aspirating components is minimized. Secure all components used in the patient’s mouth against aspiration and swallowing.
- Comply with the indicated torque levels (24 Ncm).

Adverse effects

None known.

Intended Use

Please make sure to read our comprehensive user manual for the Cercon® art software and specifically the chapter on Cercon® art custom abutments.

A) Designing and customizing the abutment

The **XIVE® custom abutment Ti by Compartis®** is to be used as delivered.

The **XIVE® custom abutment ZrO₂ by Compartis®** is to be used as delivered. The design created by the Cercon® art software already reflects the ideal anatomical shape. The **XIVE® custom abutment ZrO₂ by Compartis®** must not be reduced in order to ensure the minimum wall thickness and the minimum cementing surface are maintained.

Minimum height 4 mm (between the “preparation margin” of the abutment and it is incisal edge)

Maximum height See contraindications. If the height exceeds this maximum, the lever forces acting on the abutment would become too large, causing abutment overloading.

Screw access canal The housing of the screw-access canal is a determining factor for the longevity of the abutment. The dimensions provided by the design software must never be reduced by subsequent grinding.

B) Fabricating the superstructure

Framework design (for Cercon® base zirconia frameworks)

Minimum wall thickness 0-4 mm

Shape The design of the framework should reflect a reduced anatomic shape to support the oblique uniform thickness of the veneering ceramic layer.

Cementing

Minimum requirements for bonding areas

There is a general rule: The greater the adhesive surface, the greater the bonding area. For this reason, the abutment should be designed for its maximum allowable geometric shape, making sure, however, that the zirconia framework (crown or abutment coping) has a minimum wall thickness of 0,4 mm.

Bonding gap dimensions The width of the bonding gap must not exceed 50 µm.

Preparing the bonding surfaces All adhesive surfaces should be air-abraded with alumina (grain size, 50 µm; pressure, 2 bar).

Adhesive surface properties All adhesive surfaces must be clean and free of any greasy residue. Steam-cleaning and cleaning with alcohol is mandatory.

Cementing recommendations

The restorations can be cemented with dual-curing or self-curing composite resin cements, which should preferably be translucent in the case of all-ceramic restorations. Thoroughly remove excess cement at the crown margin.

- SmartCem 2 (Dentsply, Konstanz, Germany)
- Panavia F (Kuraray, Osaka, Japan)
- Panavia 21 TC (Kuraray, Osaka, Japan)

Miscellaneous

- The sulcus area of the abutment can be individually polished using rotating instruments. Please use the lab analogue to avoid any damage to the abutment’s implant connector.
- The abutment should preferably be restored using veneered superstructures based on our Cercon® smart ceramics zirconia system.
- Follow the recommendations given by the manufacture of the all – ceramic restorations sys-tem.

Placement

(Use original tools only)

Unscrow and remove the healing abutment. Clean the inner connection with an ear/water aerosol and dry. Connect the abutment at a torque of 24 Ncm. To tighten the screw, use the prosthetic ratchet with a 1.22 mm hex torque control device or a torque-controlled angled handpiece with 1.22 mm hex screwdriver bits. Check the fit of the superstructure on the abutment. Clean and disinfect the superstructure and cement in place.

Presentation – Storage – Disinfection

The **XIVE® custom abutment ZrO₂ by Compartis®** is not available in all countries.

The **XIVE® custom abutment ZrO₂ by Compartis®** is delivered non-sterile and is intended for single use only. Store in the original packaging until used.

CAUTION: Components delivered non-sterile must be cleaned and disinfected before they are used on patients. We recommend the following procedure:

- Cleaning:** Place the components in a cleaning bath* and clean thoroughly with a soft brush.
- Disinfection:** Soak the components in a disinfection bath* then, rinse with a. Aqua purificata and dry with a lint-free single-use tissue.
 - Follow the instructions given by the manufacture of the instrument cleaner/disinfectant with regard to concentrations and soaking time, as well as temperatures where appropriate. Only neutral disinfectants containing no chlorine, ammonia or aldehydes and with proven efficacy against HBV, HCV and HIV may be used. These must comply with all applicable provisions for disinfectants. Only freshly produced solutions may be used.
- Rinsing and drying:** Rinse components three times with distilled water (low germ count), e.g. Aqua purificata. Dry thoroughly with lint-free single-use tissues. Use oil-free compressed air to dry any bores.

4. Storage: Store the disinfected components at room temperature in a dry and dust-free location.

Copyright and trademarks

All rights reserved. No part of these Instructions for Use may be reproduced or processed, duplicated or otherwise distributed using electronic data storage and retrieval systems in any way without prior written consent on the part of the manufacturer.

Cercon®, Cercon® art, Compartis® and XIVE® are registered trademarks of DequDent GmbH, Hanau, Germany or of FRIADENT GmbH, Mannheim, Germany.

Clinical safety

The above recommendations are based on comprehensive studies including both static and dynamic (thermal cycling) loads. Please make sure to follow these recommendations to achieve a maximum of clinical safety.

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